

## 3400 SERIES

### Floor Model Universal Testing Machines

The 3400 Series universal testing machines are designed to perform tensile, compression, flex, peel, puncture, friction, shear tests, and more. The systems are compatible with hundreds of grips and fixtures found in Instron's expansive accessories catalog, with specific configurations designed to perform many of the most popular ASTM and ISO tests.

#### FEATURES AND BENEFITS

For higher force applications, the 3400 floor model systems provides up to 300 kN capacity available in short base and tall base options. Powered by Bluehill® Universal, the 3400 Series systems are simple, smarter, and safer than ever before, featuring:

- **Built to Perform:** The frame has been fortified against shock and vibration which allows for continuous testing at the maximum rated frame capacity and the maintenance free brushless AC servomotors support continuous cyclic, creep, and relaxation testing for up to 10 days.
- **Built to Last:** Reduction in debris ingress has been achieved thru gasketing and a patent-pending airflow design. Internal electrical and mechanical hardware is located safely away from directly beneath the test area. The 12 mm metal work surface is treated with an abrasion-resistant coating for durability.
- **Ergonomic Design:** The tall base option provides an ergonomically comfortable testing station. A contoured design allows the operator easy access to the test space for fixture assembly as well as a comfortable reach to the control panel and touchscreen.
- **User Safety Enhancements:** Patent pending Operator Protect architecture keeps operators safer by controlling system status and speeds throughout setup and testing via intelligent workflows. Safety Coaching visually alerts the operator to when the machine is in a safe setup mode and clearly identifies when the safety limits have been removed for testing.
- **Collision Mitigation:** Reduce damage to equipment and delicate specimens by stopping crosshead motion if force is detected upon jog or return.
- **Unmatched Support:** We service and calibrate more than 30,000 Instron systems in active use worldwide every year.



SPECIFICATIONS

		34FM-100	34FM-300
Force Capacity	kN	100	300
	lbf	22480	67440
Crosshead Travel	mm	1407	1359
	in	55.3	53.5
Vertical Test Space (A)	mm	1494	1446
	in	58.8	56.9
Horizontal Test Space (B)	mm	575	575
	in	22.6	22.6
Maximum Speed	mm/min	508	508
	in/min	20	20
Minimum Speed	mm/min	0.00005	0.00005
	in/min	0.000002	0.000002
Maximum Return Speed	mm/min	600	508
	in/min	23.6	20.0
Footprint Dimensions (h x w x d)*	mm	2287 (B1), 2587 (B2) x 1140 x 786	
	in	90 (B1), 101.85 (B2) x 44.9 x 30.9	
Position Control Resolution	nm	2.215625	1.140395
	µin	0.087	0.044
Frame Axial Stiffness	kN/mm	300	350
	lb/in	1713044	2000000
Maximum Force at Full Speed	kN	50	125
	lbf	11240	28100
Maximum Speed at Full Force	mm/min	254	153
	in/min	10	6
Weight	kg	786.2 (B1), 848.4 (B2)	858.2 (B1), 920.3 (B2)
	lb	1733 (B1), 1870 (B2)	1892 (B1), 2029 (B2)
Maximum Power Requirements	VA	2400	3000

\* The footprint width is for the system only. The Operator Dashboard monitor may add up to 500 mm (20 in) to the overall width of the frame.

**Data Acquisition Rate at the PC:**

Up to 1 kHz simultaneous on force, displacement, and strain channels.

**Load Measurement Accuracy:**

± 0.5% of reading down to 1/250<sup>th</sup> of load cell capacity

± 1.0% of reading down to 1/500<sup>th</sup> of load cell capacity

**Displacement Measurement Accuracy:**

±0.02 mm or 0.15% of displacement (whichever is greater)

**Testing Speed Accuracy:**

(Zero or constant load) ±0.2% of set speed.

**Single Phase Voltage:**

208 (-5%) to 240 (+10%) VAC, 47 to 63 Hz.

**Three Phase Voltage :**

208 to 240 VAC ±10%, 47 to 63 Hz. Compatible with both 3-phase, 4-wire, wye and delta configurations

**Operating Temperature:**

+5 to +40 °C (+41 to +104 °F)

**Storage Temperature:**

-25 to +55 °C (-13 to +131 °F)



**Humidity Range:**

+10 to +90%, non-condensing at 20 °C

**Ingress Protection (IP) Rating:**

IP 2X. Protective measures may be required if excessive dust, corrosive fumes, electromagnetic fields, or hazardous conditions are encountered.

These specifications were developed in accordance with Instron's standard procedures and are subject to change without notice. All systems conform to all relevant European standards and carry a CE mark.

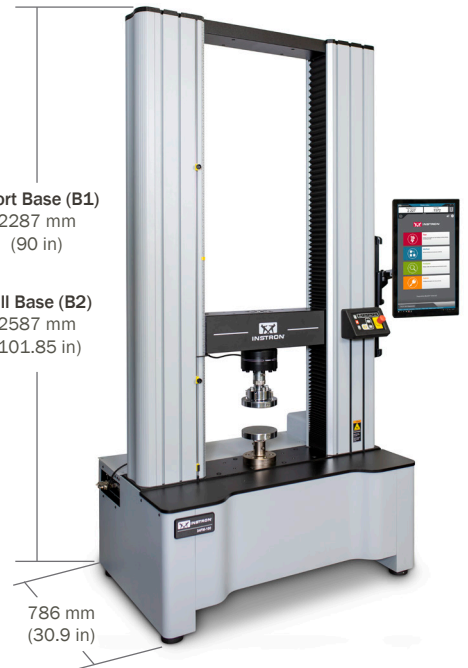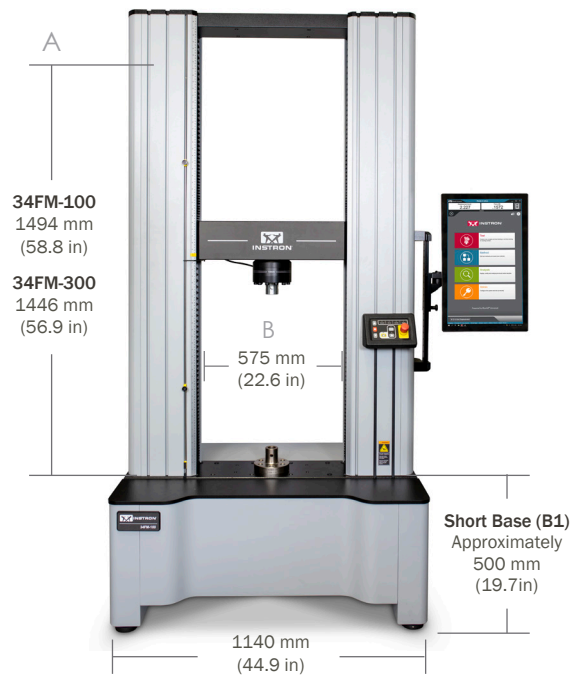

www.instron.com



Worldwide Headquarters  
825 University Ave, Norwood, MA 02062-2643, USA  
Tel: +1 800 564 8378 or +1 781 575 5000

European Headquarters  
Coronation Road, High Wycombe, Bucks HP12 3SY, UK  
Tel: +44 1494 464646